# 3D matrix for mesenchymal stem cells (MSC)

Autologous graft of MSC is used for the treatment of cardiac ischemia. The organ functionnal recovery is mainly attributed to paracrine effects, thus depends on cell survival rate after introduction. Current injection methods offer low survival rate after 3 days, limiting the efficacy of the therapeutic strategy.



## □ COMPETITIVE ADVANTAGES

- Compatible with in vivo use (biocompatible, biodegradable, ingredients and GMP process compatible)
- Adaptable and reproducible functionnality (control of 3D architecture and porosity)
- Easy handling (mechanical resistance)

## **☑** DESCRIPTION\*


- 3D matrix made of alginate and chitosan complexes suitable for MSC culture and displaying the following in vitro features:
  - Preservation of cell phenotype
  - Preservation of MSC viability to c.a. 80% after 1 month
  - Preservation, or even increase, of the paracrine secretory capacity of the MSC compared to 2D cultures
- Production process compatible to GMP standards while controlling the architecture
- property: excellent biocompatibility, significant improvements of the cardiac function and of the neovascularization of infarcted limitation the tissue, of the ventricular fibrosis

### **APPLICATIONS**

- · In vitro cell culture
- Cell therapy based on paracrine effect, especially in the field of cardiology

# () INTELLECTUAL PROPERTY


· Patent pending

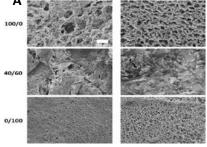

Photos : CIRIMAT & I2MC.

Matrix covering
the ischemia area

Ligation of the left
coronary artery

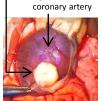

A) Surface and transversal section of the matrix porosity under various alginate/chitosan mass ratio observed by scanning electron miscroscope;

B) Picture of a matrix positionned on a rat heart 4 weeks after implantation.

# **O DEVELOPMENT STAGE**

Technology validated in lab



#### Q LABORATORIES

**CONTACT** 

- Inter-university Material Research and Engineering Centre (CIRIMAT)
- Institute of Cardiovascular and Metabolic Diseases (I2MC)





#### **ETECHNICAL SPECIFICATIONS**

Pore size	0,2μm to 400μm
Porous volume	60% to 98% of total volume
Elasticity module at 50% of deformation	1 to 100kPa

<sup>\*</sup>Technology requiring license rights.

T. +33 (0)5 62 25 50 60 sante@toulouse-tech-transfer.com www.toulouse-tech-transfer.com